



Front Panel

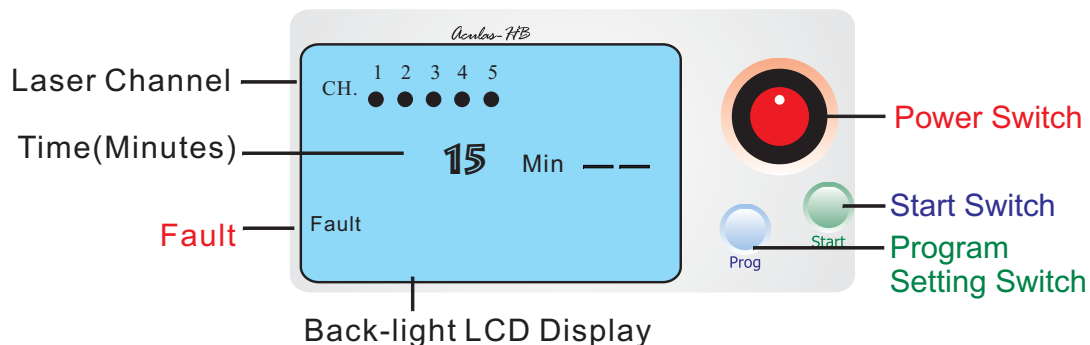

Rear Panel

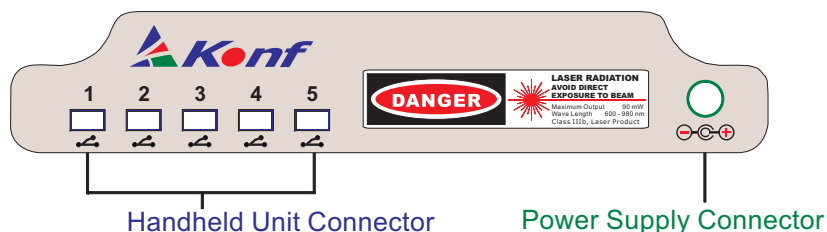

Operation			LCD Display
STEP	ACTION	PUSH	
1	Power ON		
2	Program Setting		
3	Start		
Fault Alarm:(1) Fault Sign Flash (2) Faulty Channel(s) Flash Check Connection Wire of Flashing Channel(s)			

1. Never look directly into laser beams
2. 15 minutes every 2 days is best
3. Hold on one area 5-10 seconds then move
4. Concentrate on problem areas, but expose entire area
5. There is no benefit of over exposure or more frequent sessions.

To start, turn on Power Button, wait a moment, press P button to load Program, wait for 5 circles to appear and then press S to Start.