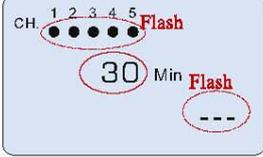
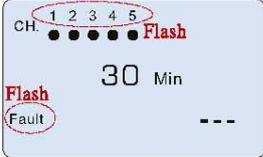


Operation

Operation			LCD Display
STEP	ACTION	PUSH	
1	Power ON		
2	Program Setting		
3	Start		
Fault Alarm:(1) Fault Sign Flash (2) Faulty Channel(s) Flash Check Connection Wire of Flashing Channel(s)			

NOTES ON USE:

1. DO NOT LOOK DIRECTLY INTO LASER, DO NOT POINT LASERS AT FACE OR EYES. LOOKING INTO THE LASER BEAMS WILL DAMAGE YOUR EYES. THERE IS NO NEED FOR SAFETY GLASSES DURING TREATMENT, BUT IF YOU DO WANT TO USE EYE PROTECTION, NORMAL UV400 SUNGLASSES ARE ADEQUATE FOR OCCASIONAL PROTECTION AS THESE WILL ALSO REDUCE VISIBLE LASER LIGHT.
2. DO NOT DISCONNECT HAND DEVICE WHILE CONTROLLER IS TURN ON, THIS MAY DAMAGE THE LASER DIODES BEYOND REPAIR AND IS NOT COVERED BY WARRANTY.
3. **PROGRAMMING IS PRESET AT 15 MINUTES FOR THIS MODEL. JUST PRESS ON, "P" BUTTON TO LOAD THE PRESET PROGRAM AND THEN PRESS "S" BUTTON TO START THE LASERS**
4. TREATMENT EVERY SECOND DAY WILL GIVE BEST RESULTS IN AROUND 3 MONTHS TIME. CONTINUAL USE IS REQUIRED TO MAINTAIN HAIR REGROWTH LEVELS.
5. IT IS POSSIBLE TO USE TREATMENT FOR MAXIMUM OF 30 MINUTES EVERY DAY, BUT 15 MINUTES EVERY SECOND DAY SEEMS BEST. ALTHOUGH IT IS NOT POSSIBLE TO BURN YOURSELF WITH LLLT LASER, IT MAY BE POSSIBLE TO OVER EXPOSE AND OVER STIMULATE HAIR FOLLICLES, SO DO NOT USE LONGER THAN RECOMMENDED OR MAXIMUM RECOMMENDATIONS ABOVE.
6. IT IS OK TO HOLD LASER HEAD DIRECTLY ON YOUR SKIN, YOU CAN EITHER MOVE THE LASER HEAD OVER THE AFFECTED AREAS SLOWLY, OR YOU CAN HOLD IN ONE AREA FOR SEVERAL MINUTES AND THEN MOVE TO NEXT AREA. TOTAL TREATMENT OF YOUR SCALP SHOULD BE 15 OR 30 MINUTES. TO USE 30 MINUTES PLEASE USE 2 TREATMENT SESSIONS IN A ROW.
7. THERE IS NO NEED TO COMB OR BRUSH YOUR HAIR. LASERS WITH LESSER POWER OFTEN USE COMBS AS THEY ARE NOT VERY POWERFUL, ACULAS AND SALON LASERS DO NOT REQUIRE COMBS TO PART YOUR HAIR.

It is normal for users of hair loss laser combs to experience an increase in hair loss for a short period after starting treatment usually around 1-2 weeks. This is thought to be caused by the stimulation of hair follicles increasing their growth rate and dropping hairs that would have otherwise been dropped soon anyway.

If excessive hair loss is developed and persists as a result of any treatment, discontinue use immediately.